



FORMALDEHYDE DISPENSING STATION

INTRODUCTION

The TDI EUROAIRE Formaldehyde dispensing Station is a compact work cabin designed for the storage, containment and dispensing of chemicals, in particular formaldehyde, in operating theatres or pathology laboratories. It is designed to be customized in dimensions and endowment for each type of requirement.

It is designed to give total protection to the operator.

WORKING PRINCIPLE OF THE STATION

- STATION – FORM is designed around our FAGV gas extraction materials, ensuring maximum protection for the user.
- By design it facilitates full access to the work area.
- Equipped with high performance motors and adjustable speed.
- Activated carbon filter (up to 10 Kg) specific for FORMALDEHIDO with high holding capacity.
- The customer's original chemical container (up to 25L) is placed on the bottom shelf and connected to a pumping system that by pressing a switch button or pedal, pumps formaldehyde into the work area faucet.



FORMALDEHYDE DISPENSING STATION





MAIN OFFICE

José Echegaray, 8 Local 1-3
28100 – Alcobendas – Madrid - Spain
Tel. +34 916 624 473 Fax.+34 916 624 469
www.tdi.es ventas@tdi.es

DELEGATIONS

Catalunya. Tel. +34 676 351 264
Galicia. Tel. +34 639 486 955
Andalucía. Tel. +34 639 481 756

FORMALDEHYDE DISPENSING STATION

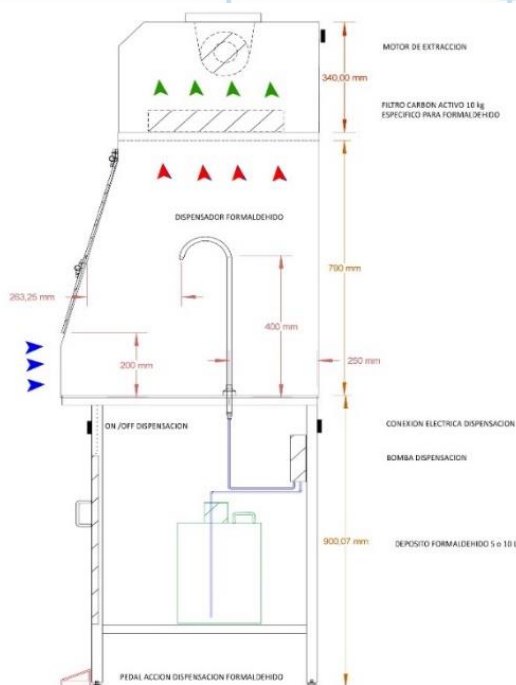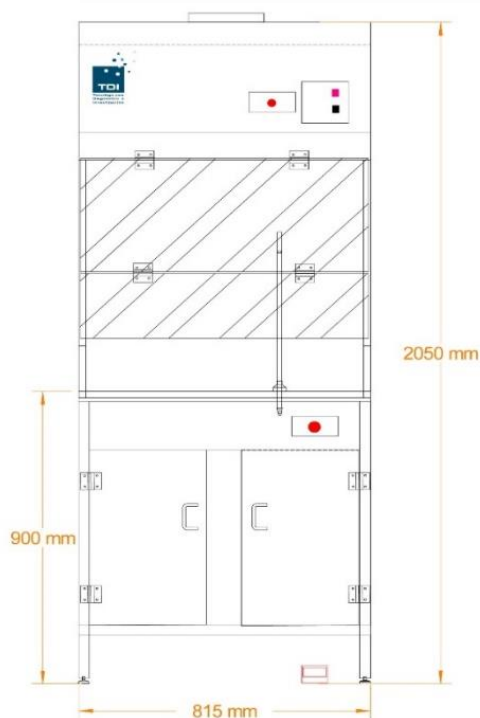

MODEL	00-ESTACION-FORM
Work area (mm) (High x Width x Depth)	790 X 800 X 25
External dimensions (mm) (High x Width x Depth)	1050 X 815 X 730
Height Surface Work	900 mm
Maximum opening two steps	400 mm y 770 mm
Barrier velocity	0,3 ---- 0,8 m/s
Noise	</= 60 dB
Outlet	PVC maximum lenght 4 meters
Outlet diamenter OPTIONAL	250 mm
Net Weight (kg)	180
Material	Exterior of epoxy painted steel sheets. Surface Stainless Steel AISI 304
Control	Aspiration speed regulator
Filter	CARBON ACTIVATED FOR FORMALDEHYDE
Fluorecent lamp	LED
Powe consumption	400 W
Powe supply	220 v, 50 h z± 10% ; 110 v , 60 Hz) Protection IP44
Formaldehido Pumping	Pump for formolin with flow max 1.3l/hour. Adjustable
Deposit formaldehde	Capacity up to 25 liters. Includes stopper with aeration valve.
Pumping Action	Includes Pedal to operate the pumping.

ORDERING REFERENCE

Reference	Description
00-ESTACION-FORM	Formaldehyde Dispensing Station Model EST-FOR 800

MAIN OFFICE

DELEGATIONS