

GARBAGE DISPOSAL BOOTH OF ANIMAL CAGES

CDB-ANI Model

INTRODUCTION

It is a suitable cabin for the treatment of animal litter in the SPF animal room.

It prevents contaminants generated from dumping animal garbage from harming experimenters and the work environment, providing safety protection for both.



CDB-ANI-1100

CHARACTERISTICS

1. Built with high-quality materials: the three side walls of the work area adopt an integrated structure of high-quality 304 stainless steel.
2. The cleanable internal parts adopt a rounded design, leaving no dead corners and easy to clean.
3. High tolerance: The working area has a large opening, which can accommodate all ordinary animal cages.
4. Control System: LCD touch screen design, real-time display of working status, easy to operate.
5. Ergonomic Design: Helps reduce fatigue, improves work efficiency and provides a safer working environment.

MAIN OFFICE

José Echegaray, 8 Local 1-3
28100 – Alcobendas – Madrid - Spain
Tel. +34 916 624 473 Fax.+34 916 624 469
www.tdi.es ventas@tdi.es


DELEGATIONS

Catalunya. Tel. +34 676 351 264
Galicia. Tel. +34 639 486 955
Andalucía. Tel. +34 639 481 756

GARBAGE DISPOSAL BOOTH OF ANIMAL CAGES

CDB-ANI Model

Flowchart

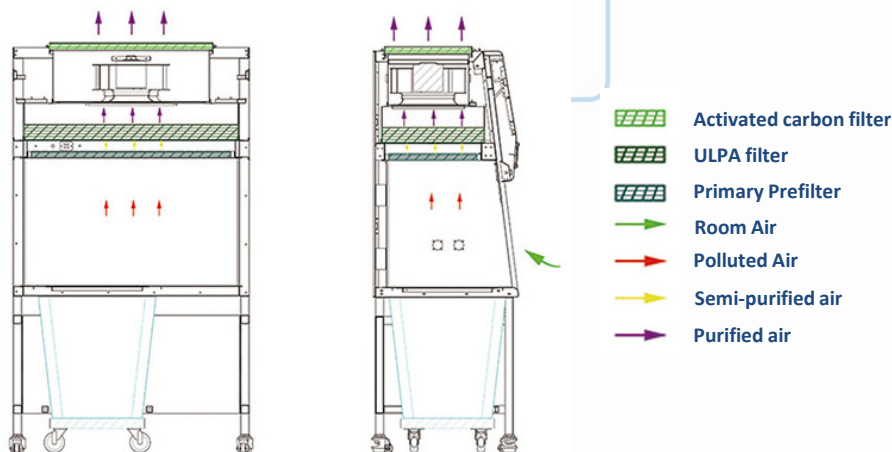


Details



Dual sensors:

The work area is equipped with temperature and humidity sensors and airflow sensors, real-time monitoring of temperature, humidity and airflow changes in the work area.

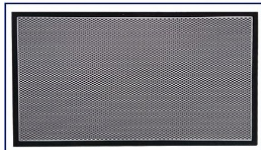

Activated carbon filter:

It has an activated carbon filter at the air outlet to absorb gaseous pollutants and reduce odor.



Large capacity storage cube:

The large capacity storage bucket is equipped with a lid and a mobile trolley, convenient for carrying.



LCD touch screen:

Shows the working time of the high-efficiency filter;
Perform the programmed time of the fan and plug;
Displays the real-time working status.

MAIN OFFICE

José Echegaray, 8 Local 1-3
28100 – Alcobendas – Madrid - Spain
Tel. +34 916 624 473 Fax. +34 916 624 469
www.tdi.es ventas@tdi.es

DELEGATIONS

Catalunya. Tel. +34 676 351 264
Galicia. Tel. +34 639 486 955
Andalucía. Tel. +34 639 481 756



GARBAGE DISPOSAL BOOTH OF ANIMAL CAGES

CDB-ANI Model

SPECIFICATIONS

MODEL	CDB-ANI-1100
Power supply	AC220v+/-, 50/60Hz (Standard); Opcional AC110V,60Hz
External dimensions (width x depth x height) (mm)	1100 X 700 X 1920
Transport dimensions (mm)	1250 X 850 X 2100
Work area dimensions (mm)	1000 X 600 X 660
Peso (kg)	240
Display	LCD TOUCH SCREEN
Laminar flow velocity	> 0,35 m/s
Noise level	< 65 dB
Lighting	> 1000 lux
LED lamps	2 X 12 W
Power	700 w
HEPA filter size	1000 X 485 X 69 mm
Activated Carbon Filter Size	758 x 408 x 35 mm
Pre filter size	967 X 453 X 30 mm
Filter efficiency	100,00%

ORDER INFORMATION

Item	Description
425-CDB-ANI-1100	Garbage Disposal booth of animal cages CDB-ANI Model



MAIN OFFICE

José Echegaray, 8 Local 1-3
28100 – Alcobendas – Madrid - Spain
Tel. +34 916 624 473 Fax.+34 916 624 469
www.tdi.es ventas@tdi.es

DELEGATIONS

Catalunya. Tel. +34 676 351 264
Galicia. Tel. +34 639 486 955
Andalucía. Tel. +34 639 481 756