


→ Animal bedding disposal cabinet - Top Flow 2



Laboratory Animal Technology: Air Flow Systems



Labor- und Pharmatechnik

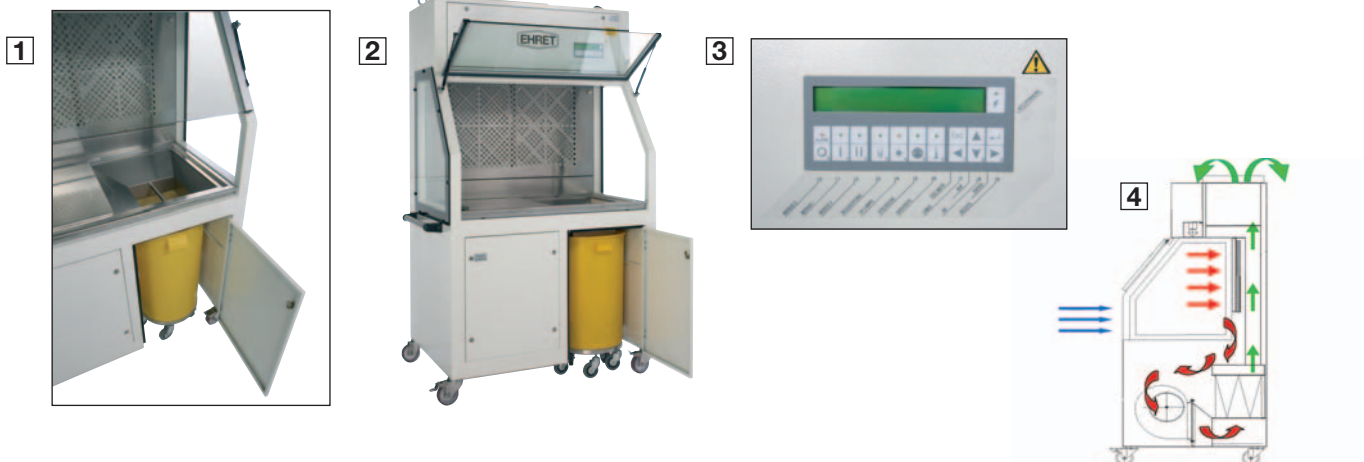

1 With waste trolley on wheels 2 With open front panel 3 Soft touch display 4 Functional layout

Top Flow 2

The Top Flow 2 is an indispensable assistant for **reduction of allergens at the work place.**

Due to the air velocity of 0,8 m/s it is granted that the particles are disposed off via the pre-filter and the HEPA filter when emptying the cages. The air which surrounds the device thereby remains free from allergens.

Manipulation and instrumentation are adapted to modern safety workbenches, adjustments and alarms are visible as plain text on the **soft touch display, which is easy to clean.**

Failures, e.g. air velocity too low, will also be signalised acoustically.

Main filter: H13 (HEPA) as per EN 1822

Exhaust air filter: H13 (HEPA) as per EN 1822 (alternatively also with active carbon filter)

Pre-filter: G3 as per EN 779

Type	Outer dimensions W x H x D	Work space W x H x D	Power	Weight	Code-no.
Top Flow 2	1270 x 2010 x 856 mm	955 x 678 x 630 mm	700 W	approx. 250 Kg	Top Flow 2